


# Emergency Power Down Procedures for the 3 Tesla Scanners

3 Tesla Core Facility




# What to Do?

- In case of power interruptions, water etc.,
  - Make sure the subject is safe
  - Make sure researchers are safe
- Notify appropriate personnel
  - NIH
  - Maintenance
  - Housekeeping
  - GE
  - 3T Core / FIM
- Power down Signa, ADW, etc (Do **NOT** Quench)

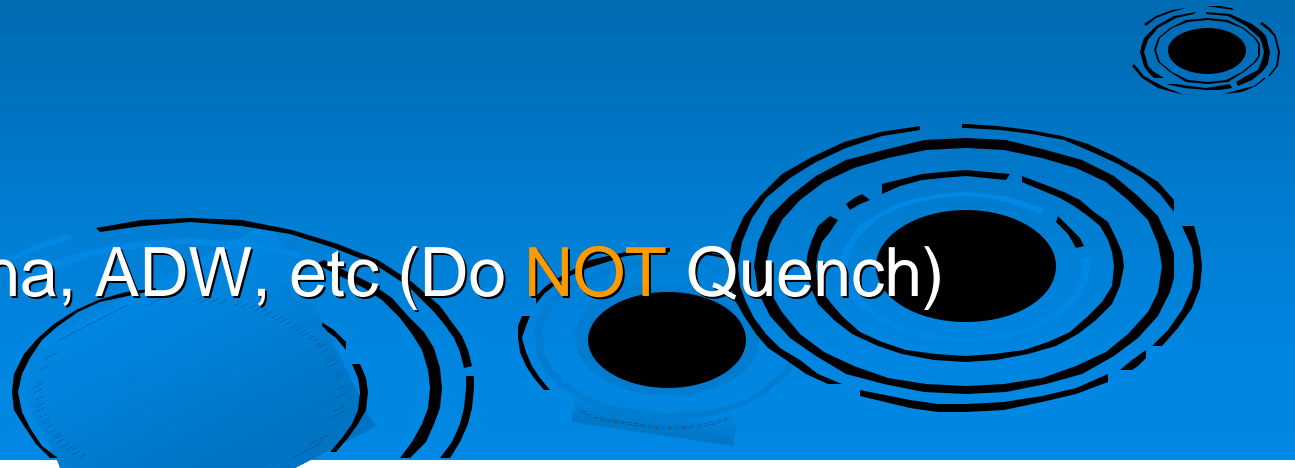


# Safe Removal of Subjects

- Remove subject from scan room prior to Signa shutdown (if possible)
- Escort the subject to the waiting area or to an alternate area of safe refuge

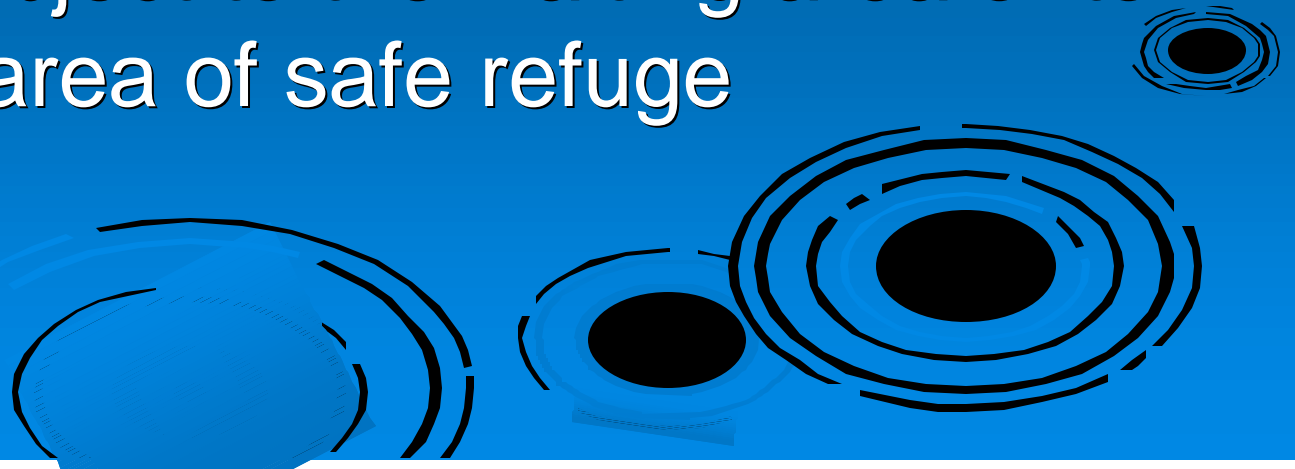


# Major Power Interruptions

- While the NIH itself is on a backup generator, it is not unknown to have long lasting power failures (2+ hours). This presentation will address power down procedures for the software, which is **NOT** the same as a magnet quench
  - NIH may need to divert power for critical patients and facility operations
  - PEPCO may need to divert power to other critical sites
- The UPS systems on Signa and peripheral computers last about 20 minutes, enough to allow for an orderly power down of software
  - Signa console
  - AW/ADW
  - Potassium / Aluminum
  - Peripheral PC's, including the power monitor and timing PC

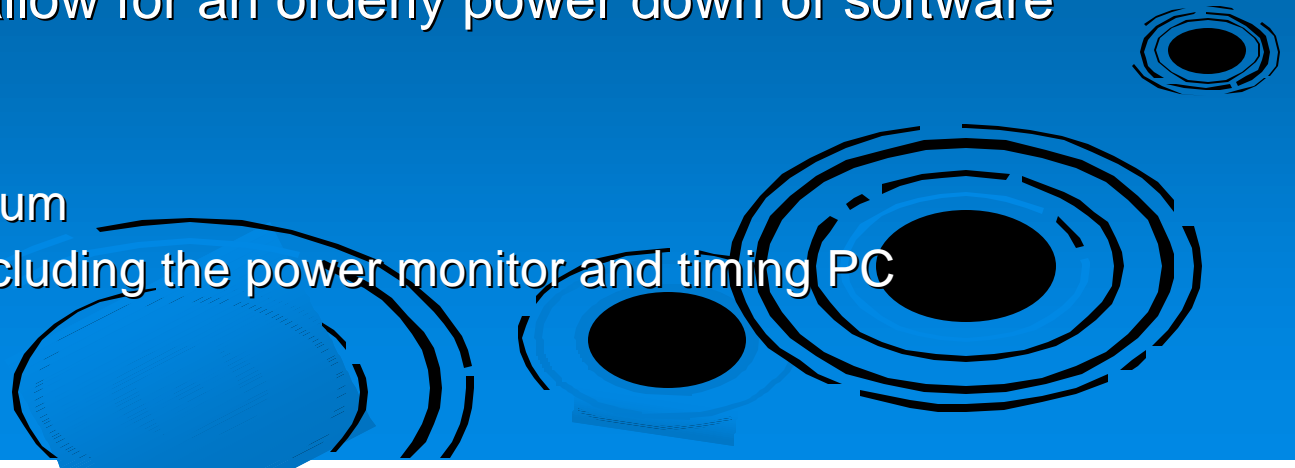

# Water from the Ceiling, Floors etc


- If you notice water in an unusual area e.g., from the ceiling, on the floors, etc....
  - Notify CC Maintenance **immediately** @ 301-496-5862
    - Be prepared to give the building number, room number, your name and a call back number
  - Call CC housekeeping @ 301-496-2417

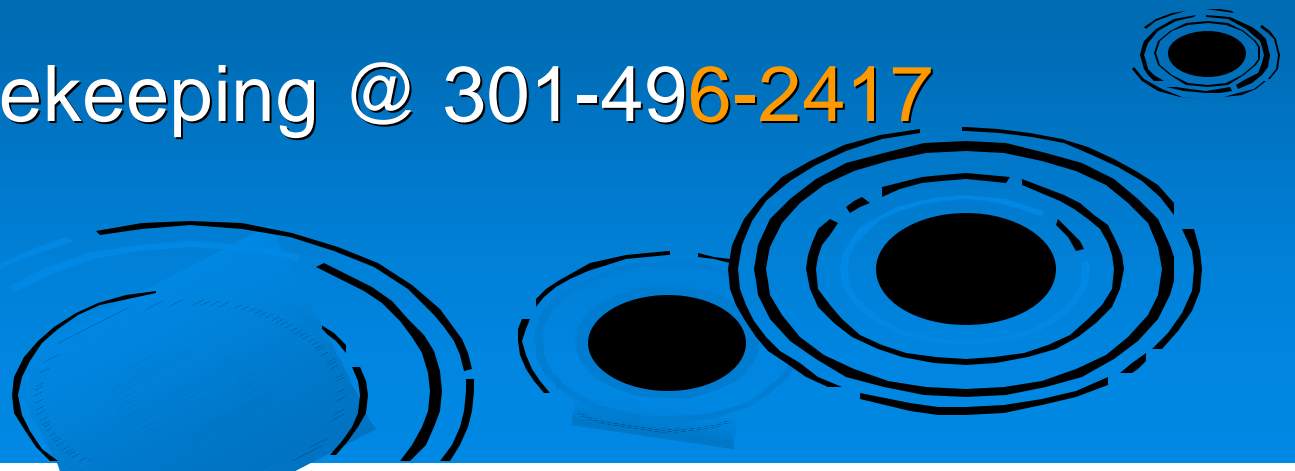

### 3 Tesla Core Facility / FIM

Name	Office	Home	Cell	Pager
Peter Bandettini	402-1333			
Sean Marrett	402-1378		301-905-8059	
Jerzy Bodurka	402-1343		301-728-0606	
Wen-Ming Luh	496-2623		301-728-0693	
Rasmus Birn	402-1350	301-519-9670	301-461-2975	
Adam Thomas	402-6351	301-915-0212	240-505-4872	
Karen Bove Bettis	435-4852	703-242-2261		104-6475
Paula Rowser	451-3877			
Others				
Brad Burge (GE)	594-4524			104-6511
Warren Embrey (GE)	594-4524			104-6527
CC Maintenance	496-5862			
CC Housekeeping	496-2417			

# Signa Shutdown End Exam

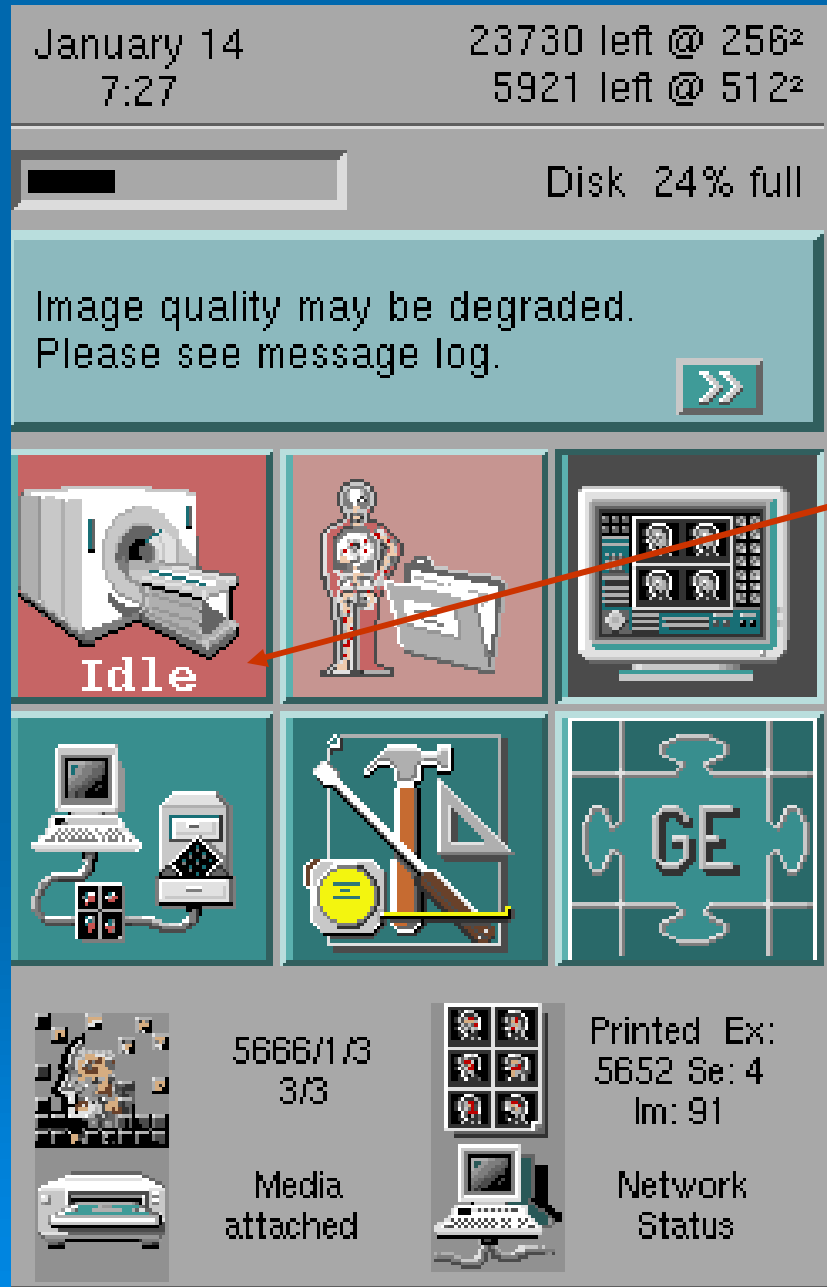

- Click on the Scanner Operations Icon (Idle)
  - Click End Exam

# Signa Shutdown


- Click on the Tools Icon
  - Click on System Shutdown
    - System Shutdown
      - Yes or No?
      - Yes





# Power Down the Sun ADW (3T1)

---

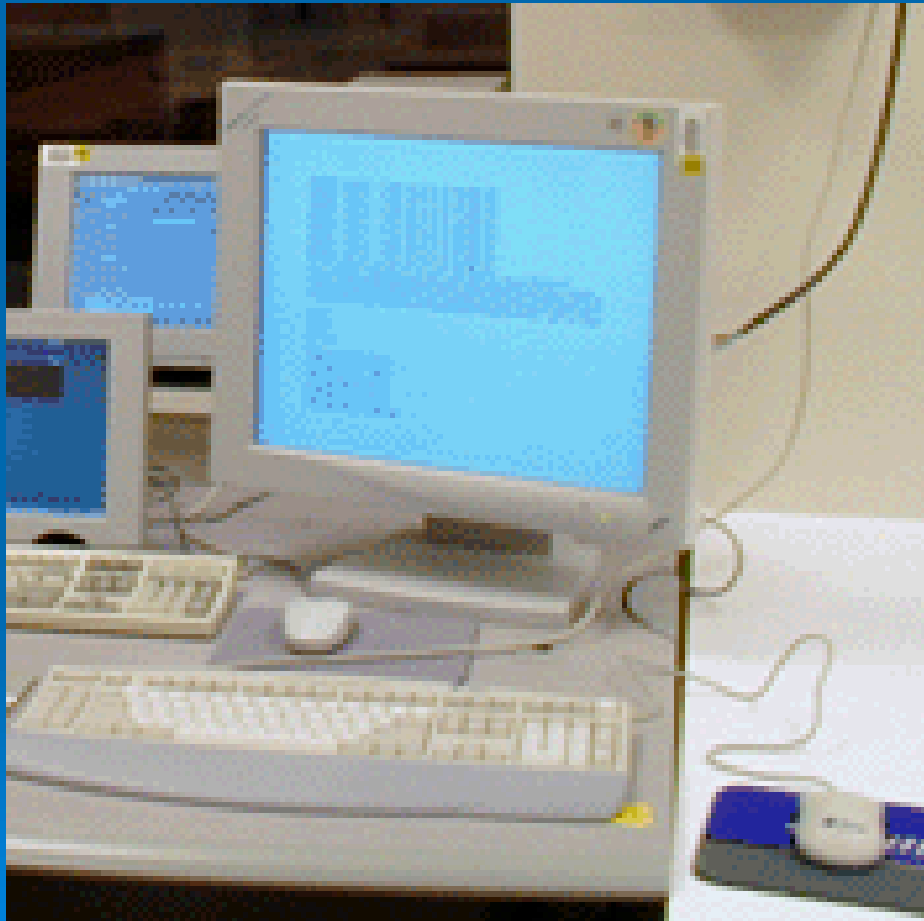

3T 1 Console

- Connected to the Signa IDBM Database
- Resides to the right of the scan console
- If there is time...save your data 1<sup>st</sup> and move to SDC Data (Potassium/Aluminum)

# SUN ADW Screen (3T1)

\* rmb = right mouse button



console@bin/csh/bin/csh

Cmdtool -/bin/csh

fim3t-adw [102]%

Using the *rmb* click on:  
Service Tools

Applications

RealTime Control  
Neuro Apps  
IDBM Apps  
Terminals  
Tools  
Window Manager  
Service Tools

# SUN ADW (3T1)

Cmdtool -/bin/csh

fm3t-adw [102]%

Using the *rmb* click on  
Logout

## Applications

RealTime Control  
Neuro Apps  
IDBM Apps  
Terminals  
Tools  
Window Manager  
Service Tools


## Service Tools

Install License Key  
Logout




console@bin/csh/bin/csh

# SUN ADW (3T1)

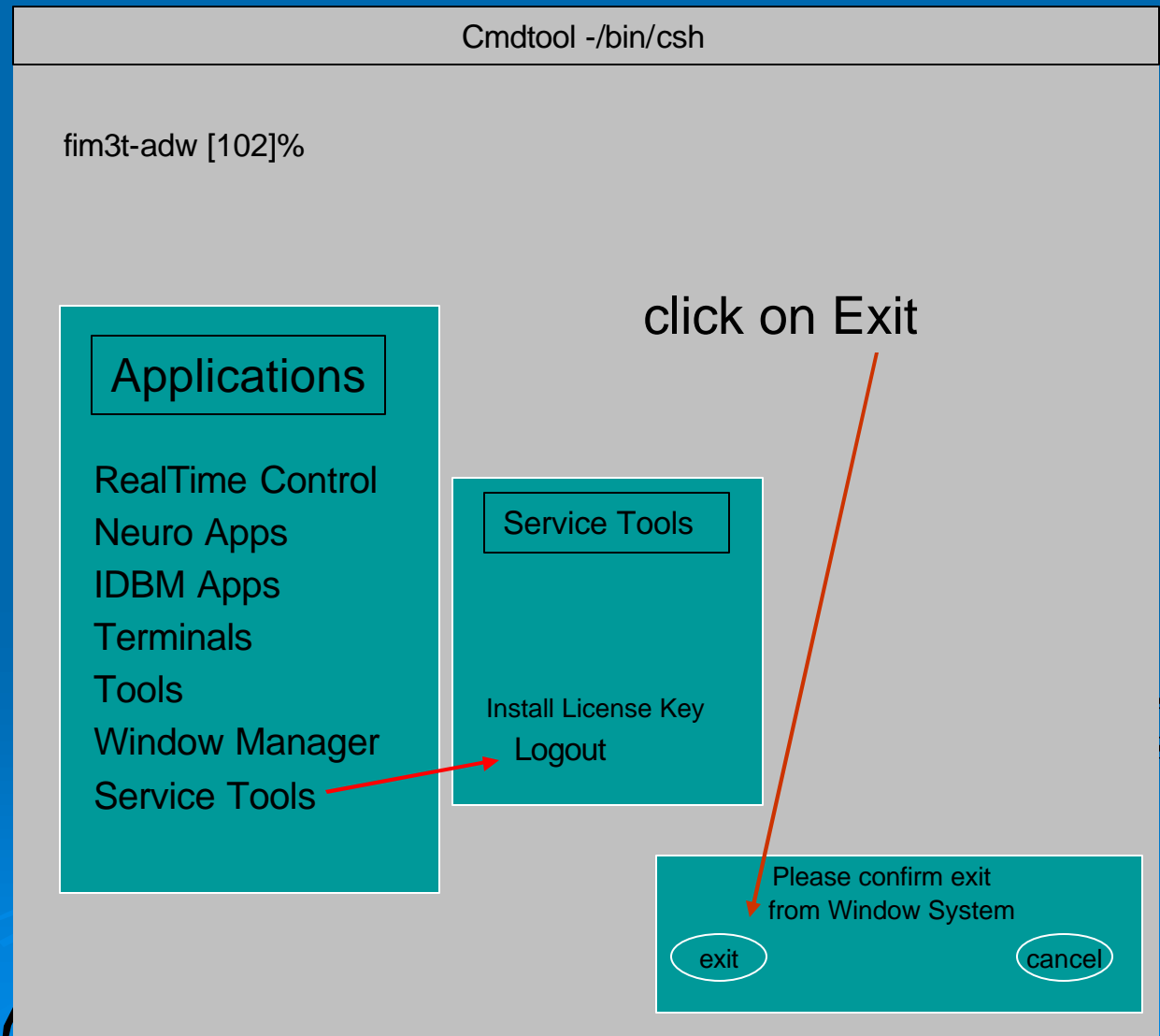

# Sun ADW (3T2)

- At the Sun/Advantage workstation in 3T2
- Click on the Security Icon
  - Click Shutdown
    - Are you sure? [yes]
      - Yes

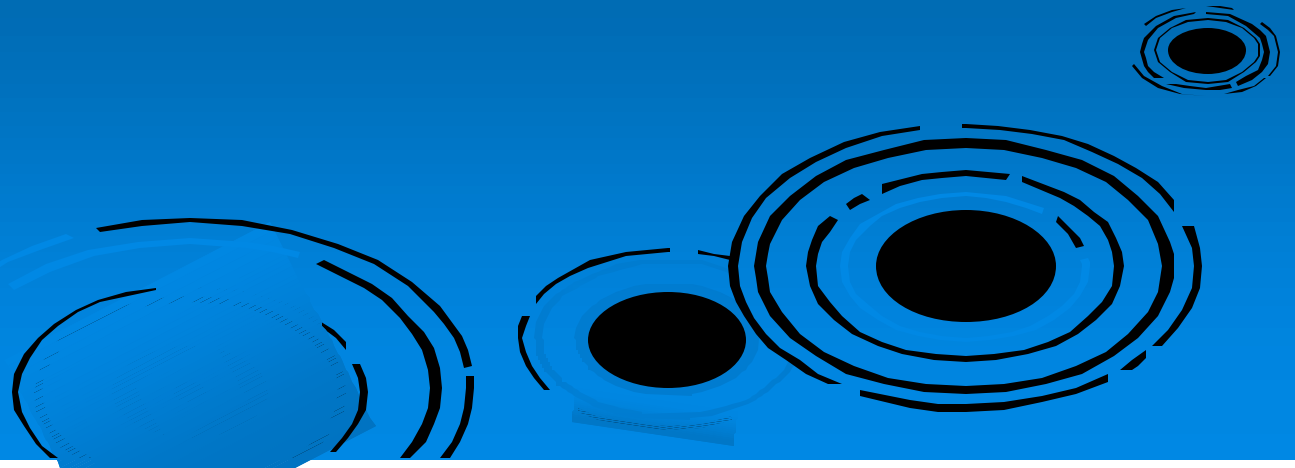

# Potassium/Aluminum Screen


(sdc-adw or sdc-nfs)

Sdc-adw@potassium.nimh.nih.gov:/home/sdc-adw-konsole

fMRI  
done

Sdc-adw@potassium.nimh.nih.gov:/home/sdc-adw-konsole

Using the *rmb* click on:  
Logout

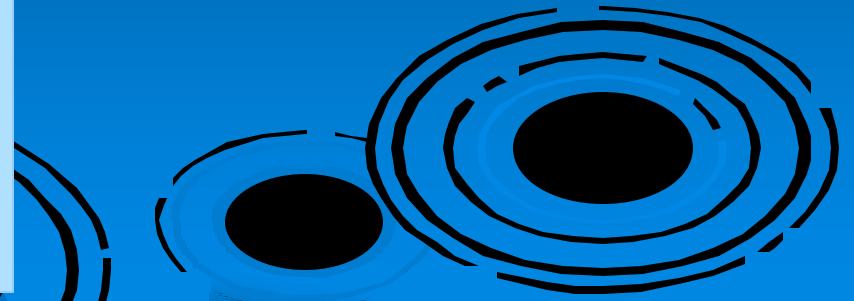

# Peripheral Computers

The peripheral computers are PC and Mac based. In 3T1, they are wall-mounted and in 3T2, they are located on the workbenches. Follow normal shutdown procedure as you would for any PC or Mac.

